

Precision-Trained Deep Learning: Redefining Cardiac Imaging

Kirsten Boedeker, PhD, DABR
Senior Manager, Medical Physics
Canon Medical Systems Corporation

Introduction

Canon Medical's new Precise IQ Engine (PIQE) algorithm delivers sharp cardiac images without compromising on noise, temporal resolution or dose. PIQE is a Super Resolution Deep Learning Reconstruction (SR-DLR) algorithm that maximizes the Aquilion ONE / PRISM Edition's 16 cm detector for extraordinary spatial resolution, reduced noise and single-rotation cardiac coverage, leading to more confident diagnosis of small vessels, plaques, and fine cardiac structures.

PIQE is trained using high quality data acquired on the commercially-available Aquilion Precision CT system, which features Ultra-High Resolution (UHR) 0.25 mm³ detectors—approximately twice the resolution of conventional detectors—in routine clinical practice.

With no loss of low contrast detectability and no additional radiation dose to the patient, PIQE optimizes spatial resolution for the most clinically relevant tasks while simultaneously decreasing noise. PIQE features a next generation, three-dimensional neural network trained to identify signal both in-plane and longitudinally throughout the cardiac volume dataset. This helps to ensure continuity of small, longitudinally running vessels, which are often obscured by conventional reconstruction algorithms.

Ultra-High Resolution CT and Deep Learning Reconstruction

Aquilion ONE / PRISM Edition

The Aquilion ONE / PRISM Edition offers One Beat Cardiac imaging: covering the whole heart in a single rotation via 320 detector rows of 0.5 mm thickness. The ability to image the heart in one 0.275 second rotation prevents misregistration and excludes artifacts caused by stitching or beat-to-beat variation and virtually eliminates contrast differences between the aortic root and apex of the heart.

PIQE builds off the foundation of Advanced intelligent Clear-IQ Engine (AiCE), a Deep Learning Reconstruction algorithm for CT. The training of PIQE harnesses the power of the Aquilion Precision Ultra-High Resolution CT scanner to enhance cardiac image quality on the Aquilion ONE / PRISM Edition.

Aquilion Precision

The Aquilion Precision CT system, introduced in 2017, was engineered from the ground up to create Ultra-High Resolution CT images. The system features 1792 detector elements per row, double the number of a conventional system, resulting in twice the intrinsic in-plane spatial resolution of a conventional CT detector. In addition, each of the 160 detector rows along the z-direction is 0.25 mm thick, half that of a standard CT detector. Proprietary cutting techniques allow for optically isolated detector elements and ultra-thin septa, resulting in a substantial increase in light-sensitive area relative to conventional CT. The Aquilion Precision detector is paired with an advanced X-ray tube design, utilizing reduced focal spot sizes, as small as 0.4 mm × 0.5 mm and rotating at 10,000 rpm to efficiently dissipate heat. Figure 1 demonstrates the spatial resolution in a cardiac case scanned with Aquilion Precision.

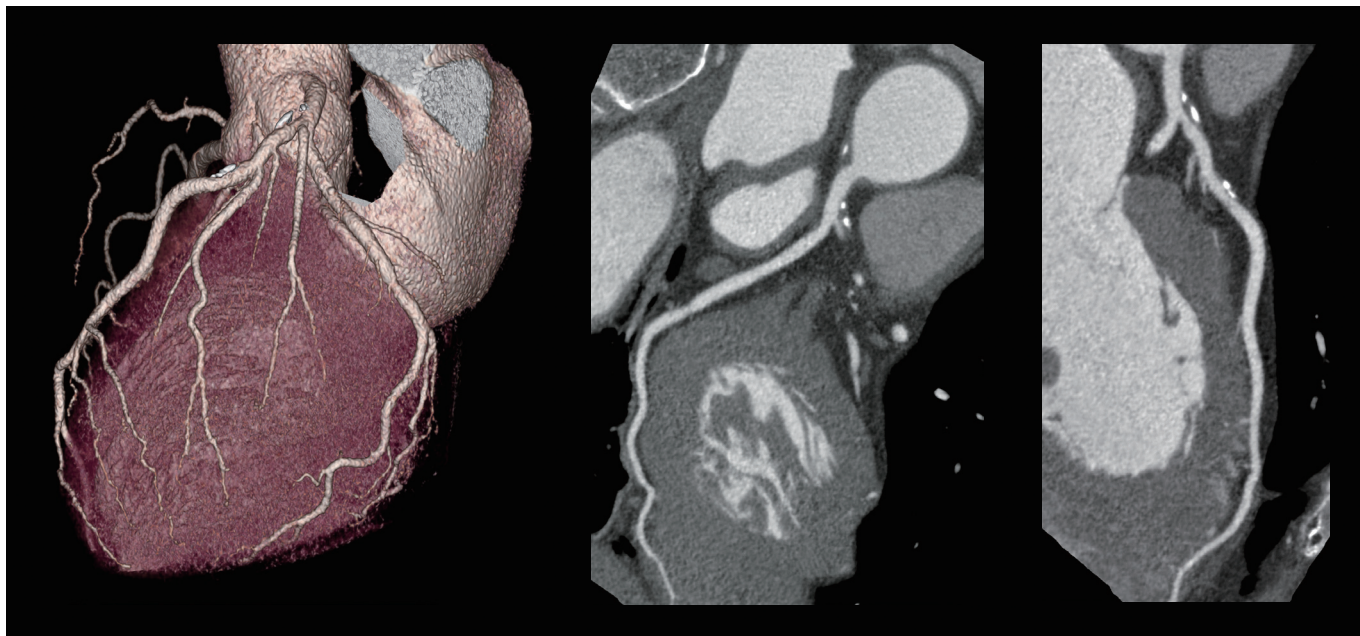


Figure 1: In this cardiac case performed on Aquilion Precision the sharp depiction of calcified and non-calcified plaques in the proximal LAD improves diagnostic confidence.

Deep Learning Reconstruction

AiCE Deep Learning Reconstruction was originally developed in order to maintain the spatial resolution benefits of the Aquilion Precision detector while achieving dose neutrality* relative to conventional resolution CT. The ability of Deep Learning to improve spatial resolution and low contrast detectability while reducing noise is due to the power of machine learning, a form of Artificial Intelligence (AI).

Deep Learning is a type of machine learning that uses multi-layered neural networks to perform a task. Neural networks are comprised of thousands of “neurons,” each of which perform a mathematical operation on the image data. The computational power behind thousands of neurons working together allows the network to generate more sophisticated rules for identifying signal features than is possible with conventionally programmed reconstruction. In the case of AiCE DLR, the task the network learns to perform is simply to distinguish signal features from noise. It identifies and enhances signal while removing noise, resulting in an industry-first low contrast specification of 1.5 mm at 3 HU and 22.6 mGy.

The key to a successful Deep Learning neural network lies in how it is trained, i.e. the process by which the neurons learn to perform the desired task. In order to learn, the network compares its output image to a gold standard reference image. The network adjusts the weights of its neurons until the error with respect to the gold standard is minimized. In the case of AiCE DLR, the gold standard clinical reference images are acquired with high tube current and reconstructed with a Model-Based Iterative Reconstruction (MBIR), that takes into account modelling of system optics, system physics, scanner statistical properties and human anatomy, and uses a greater number of iterations than could

be otherwise used in a clinical setting due to time constraints.

Aquilion Precision’s UHR mode reconstructed with AiCE DLR not only preserves more high contrast spatial resolution than traditional reconstruction, it also achieves dose neutrality relative to conventional resolution CT.

PIQE builds off the foundation of Advanced intelligent Clear-IQ Engine (AiCE), a Deep Learning Reconstruction algorithm for CT. The training of PIQE harnesses the power of the Aquilion Precision Ultra-High Resolution CT scanner to enhance cardiac image quality on the Aquilion ONE / PRISM Edition.

PIQE Overview

PIQE is trained with Aquilion Precision data in order to maximize the inherent spatial resolution on the Aquilion ONE / PRISM Edition. This type of Deep Learning neural network is called a Super Resolution algorithm. With conventional reconstruction, there has long been a tradeoff between spatial resolution and noise that has forced CT systems to operate at much lower levels of spatial resolution than the detector and focal spot can produce in order to maintain clinically acceptable dose levels. With Deep Learning this resolution can be utilized without increasing dose or noise. As with all Deep Learning, the key to maximizing performance lies in the training process.

* Applies to AiCE Body only

PIQE Training

The neural network behind PIQE is trained using high quality cardiac cases acquired on the Aquilion Precision. These cases are acquired with the Aquilion Precision's Super High Resolution (SHR) mode that yields Ultra-High Resolution in-plane and 0.25 mm nominal slice width. The SHR data is reconstructed with AiCE DLR which implicitly contains all the advanced models of MBIR. In addition to UHR mode, the Aquilion Precision also has Normal Resolution (NR) mode that combines detector channels to generate conventional resolution images in-plane and 0.5 mm nominal slice width equivalent to Aquilion ONE / PRISM Edition. Raw data acquired in UHR mode can be reconstructed through a downsampling algorithm to yield simulated Normal Resolution images that have been demonstrated to be

equivalent to true Normal Resolution images.¹ With this approach, both UHR and NR images can be produced from a single acquisition. As a result, pairs of UHR and simulated NR images have perfect spatial alignment, ideal for training a neural network. Simulated NR images are input to the neural network and the corresponding UHR images are used as the gold standard target images. The neural network learns to maximize the inherent resolution possible with NR images while decreasing noise (Figure 2).

The neural network behind PIQE does not learn features solely in the axial plane but rather in three dimensions, meaning signal features are identified and preserved in all three planes. This makes PIQE well-suited for cardiac exams, which are usually reviewed in MPR and curved planes.

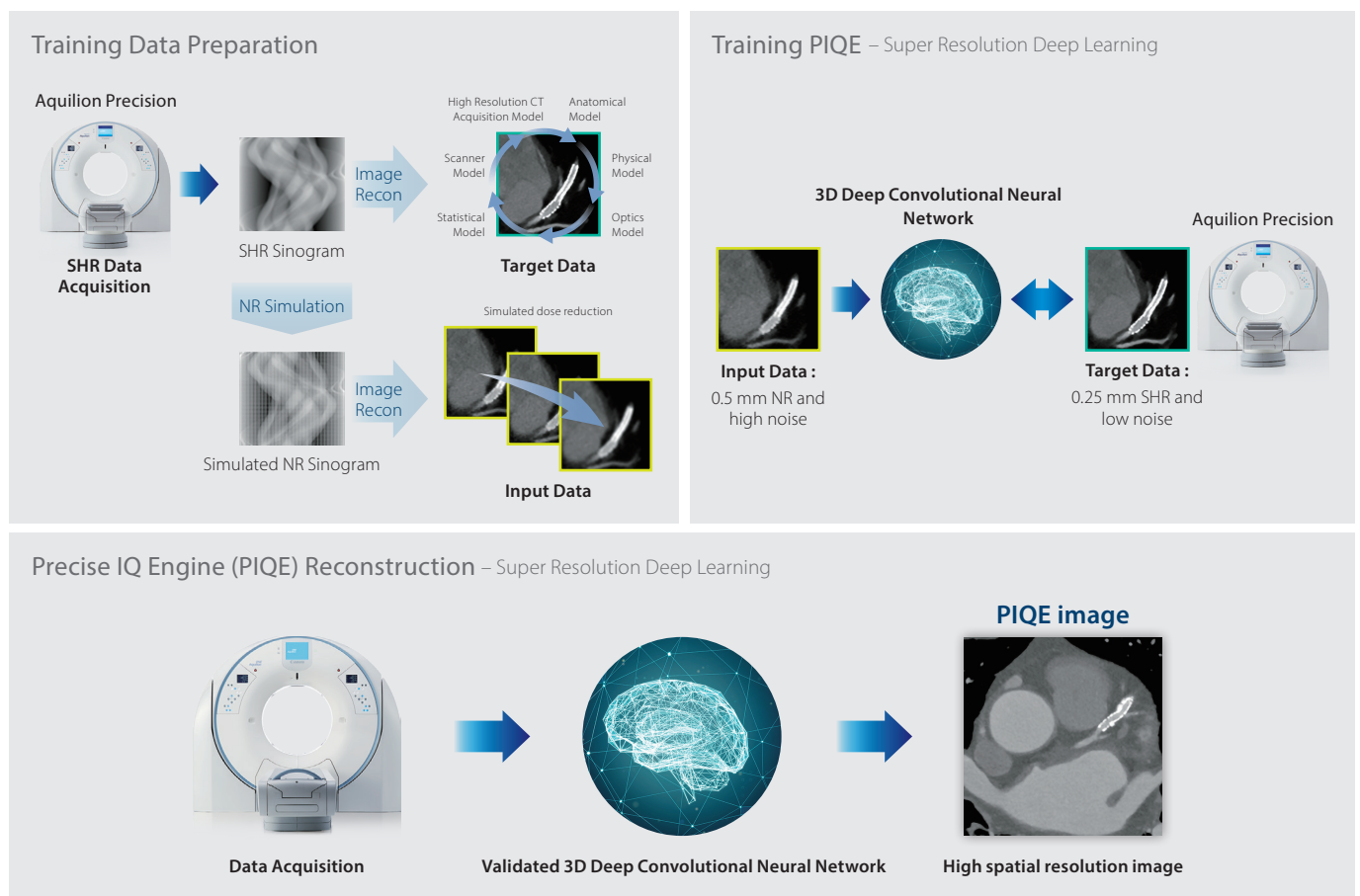


Figure 2: Training of PIQE network. Training data is prepared by down-sampling Precision UHR data. The network is trained using UHR data paired with down sampled, simulated NR data and actual UHR data. Once trained the network is validated and applied to the Aquilion ONE / PRISM Edition Scanner where it does not continue to learn.

PIQE Performance

CT is a well-established tool for the non-invasive detection and rule-out of coronary artery disease, currently the number one cause of death worldwide. PIQE increases spatial resolution by 5-7lp/cm*¹ and doubles the signal definition*² relative to conventional hybrid iterative reconstruction (AIDR*³ 3D). PIQE is optimized for routine clinical use, which means maximizing resolution in realistic, typical field of views. Furthermore, the CNR with PIQE is increased in all three dimensions improving image quality in all planes for multiplanar and curved reformats.

PIQE's improvement in spatial resolution is achieved with no loss in low contrast detectability (Figure 3).

This means objects are better defined and easier to visualize than ever before (Figure 4).

These advancements in spatial resolution foster improved visualization of small vessels, stents and in-stent structure with no dose penalty. Furthermore, the three-dimensional neural network helps promote continuity of vessels in all planes.

PIQE not only increases spatial resolution but reduces noise by 45% and yields identical low contrast compared to conventional hybrid iterative reconstruction (Figure 5 and 6). Soft, lower contrast plaque is often the most clinically relevant for patient care. Note the noise texture is not impaired by PIQE and is equivalent to conventional hybrid iterative reconstruction (Figure 5).

PIQE provides better visualization of fine structure

The spatial resolution abilities of PIQE help mitigate partial volume effects that can distort CT number accuracy on small, fine structures. PIQE demonstrates better strut delineation and HU values that are more accurate vs conventional hybrid iterative reconstruction (Figure 4).

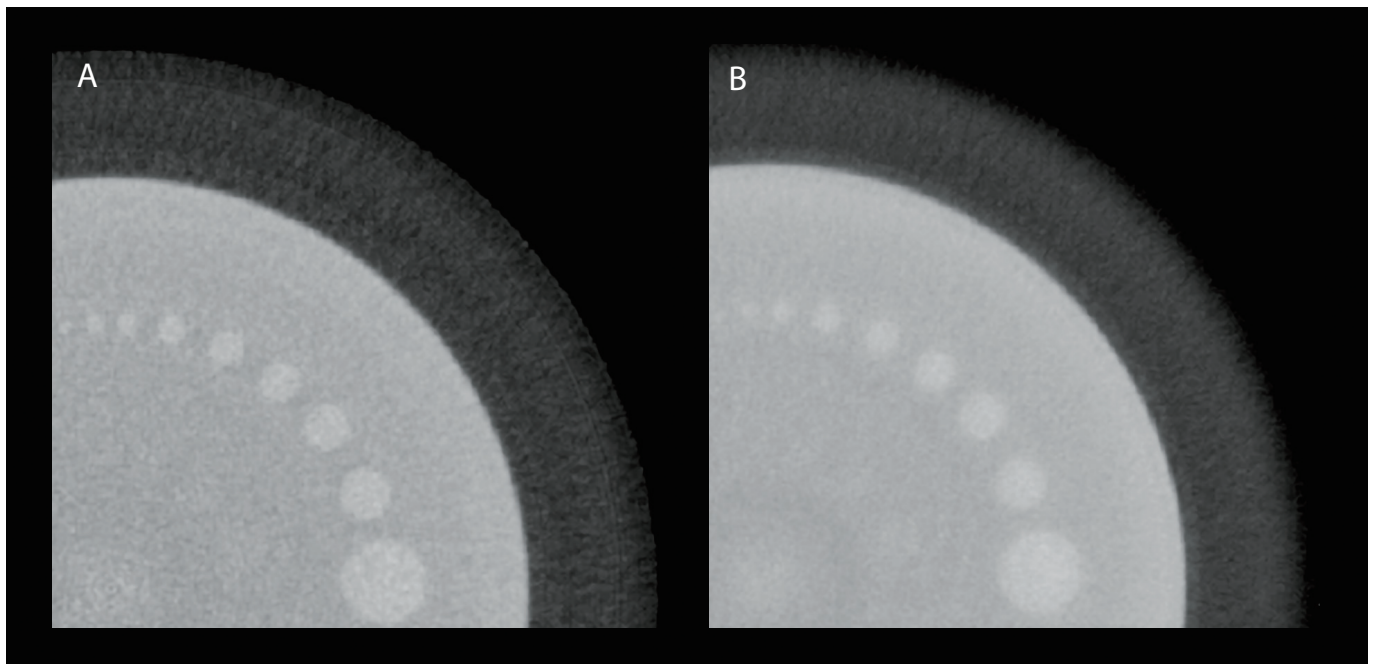


Figure 3: Catphan® phantom demonstrating equivalent low contrast detectability AIDR 3D Enhanced (A) compared with PIQE (B)

* 1 Mid to high contrast objects

* 2 Average increase in magnitude of the mid and high contrast MTF across frequency

* 3 Adaptive Iterative Dose Reduction

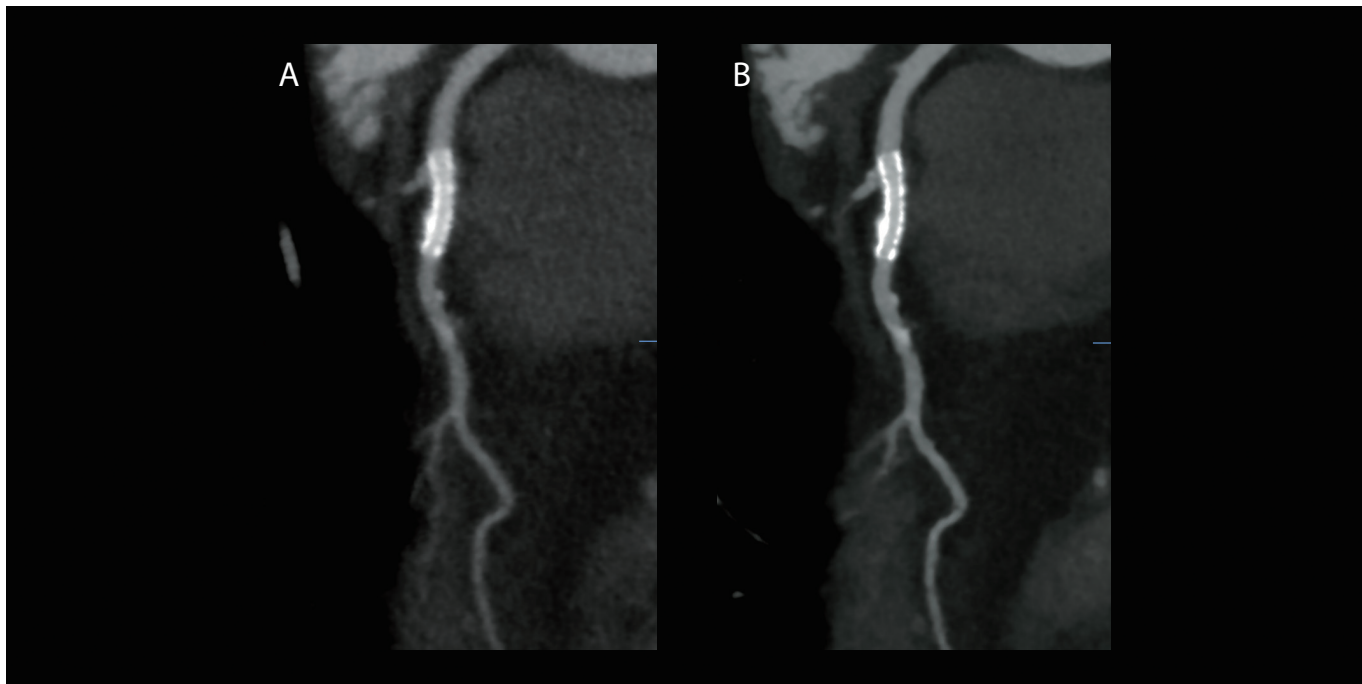


Figure 4: Visualization of a stent placed in the LAD using either AIDR 3D Enhanced (A) and PIQE (B). The stent as well as the intra-stent lumen are well defined on the PIQE images.

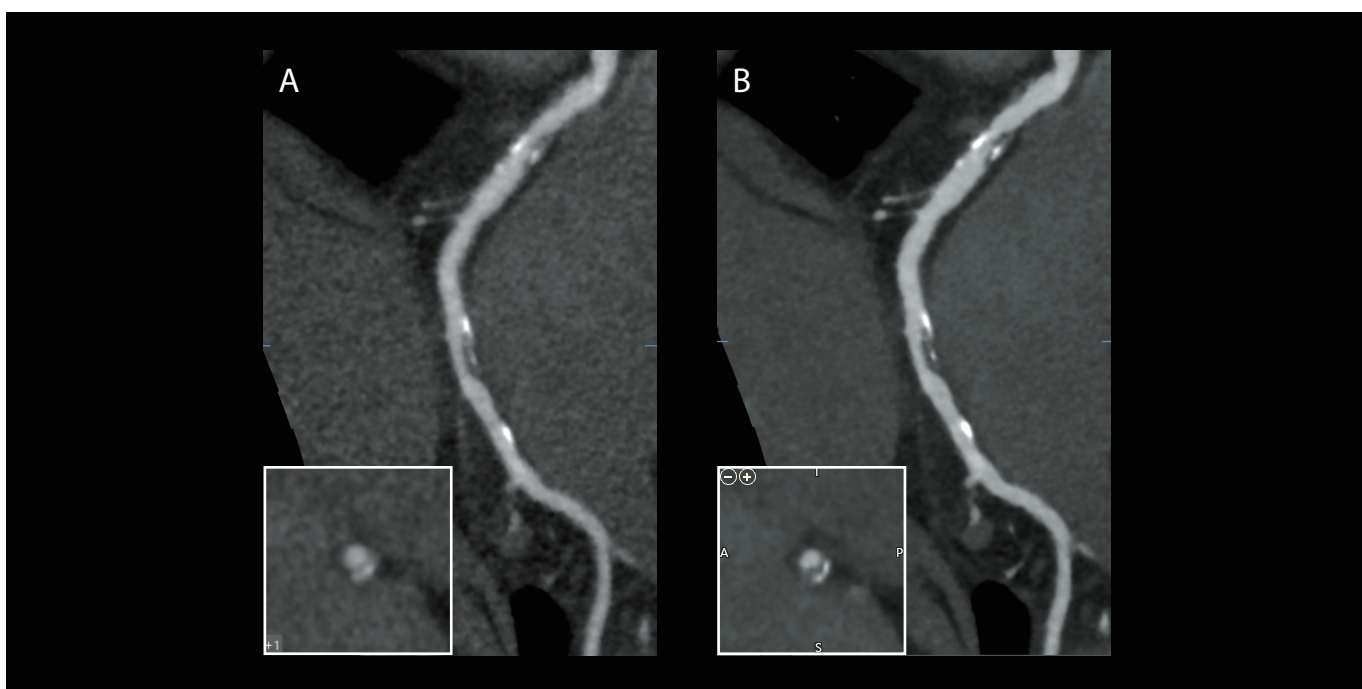


Figure 5: Multiple lesions are seen in the RCA. The various lesions and their different plaque components are clearly delineated in the PIQE images. AIDR 3D Enhanced (A), PIQE (B).

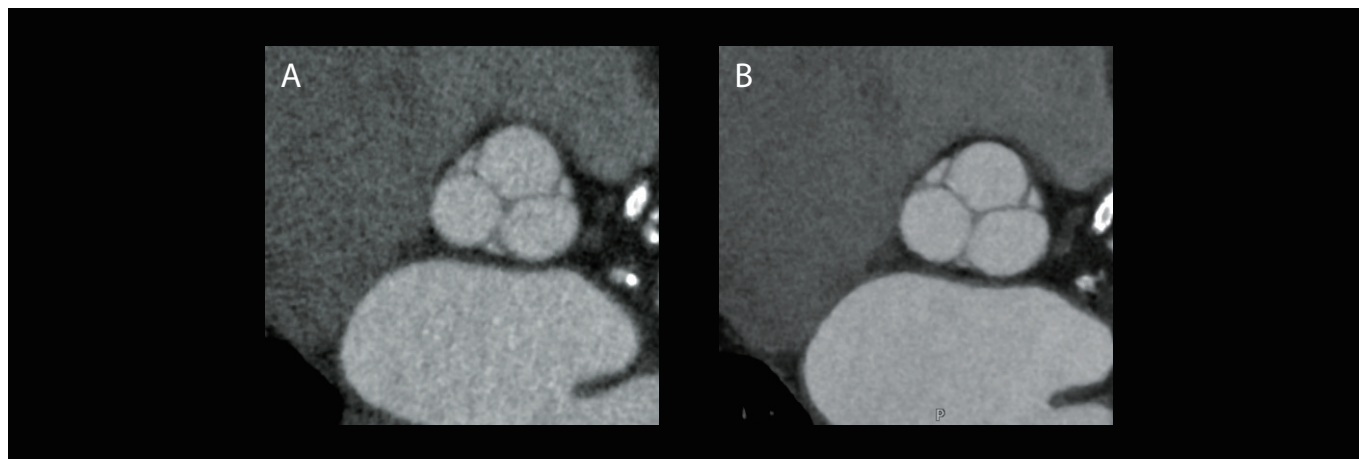


Figure 6: Aortic valve leaflets are clearly depicted in the PIQE reconstruction. AIDR 3D Enhanced (A), PIQE (B).

Cardiac Suite

PIQE joins a suite of features available on the Aquilion ONE / PRISM Edition designed to optimize cardiac image quality and workflow.

High quality cardiac scans start with patient positioning. ^{SURE}Position automatically corrects both vertical and lateral patient position and then produces a new virtual scanogram, obviating the need for a repeat scanogram, saving time and dose.

^{SURE}Cardio offers a number of tools that helps ensure consistency of image quality across patients of varying sizes and heart rates. Real Time Beat Control analyzes and predicts an individual patient's heartbeat in order to optimize ECG gated scanning. Real Time Beat Control also includes automatic Arrhythmia detection, which detects irregular heartbeats and adapts the scan protocol accordingly. PhaseXact detects the optimal cardiac phase with the least amount of motion for clinical evaluation.

^{SURE}kV and Automatic Exposure Control for automatic kVp and mA selection can be used in combination with PIQE for maximum dose efficiency.

Canon Medical's ECG-gated Single Energy Metal Artifact Reduction (SEMAR) corrects metal artifacts created by implanted cardiac devices such as pacemakers and defibrillators, to optimize the readability of cardiac CTA images. SEMAR also helps reduces artifacts due to photon starvation and scatter. SEMAR can be preset in the scan protocol such that the application of SEMAR is fully automatic during reconstruction, requiring no additional input from the operator.

With a full volume reconstruction time of under two minutes, PIQE can be integrated into the standard workflow for all cardiac examinations.

Conclusion

PIQE uses a 3D Deep Convolutional Neural Network to bring the power of Ultra-High Resolution CT to full coverage, single rotation cardiac examinations on the Aquilion ONE / PRISM Edition, with no increase in patient dose. PIQE increases spatial resolution, lowers noise, and improves CNR, all with no loss in low contrast detectability relative to conventional hybrid iterative reconstruction. The combination of improved spatial resolution while maintaining LCD helps improve visualization of stents and coronary plaque, both calcified and non-calcified, providing a more confident diagnosis of coronary artery disease.

Reference:

1. Validation of synthesized normal-resolution image data generated from high-resolution acquisitions on a commercial CT scanner.
Hernandez AM, Shin DW, Abbey CK, Seibert JA, Akino N, Goto T, Vaishnav JY, Boedeker KL, Boone JM.
Med Phys.

CANON MEDICAL SYSTEMS CORPORATION

<https://global.medical.canon>

©Canon Medical Systems Corporation 2021-2023. All rights reserved.
Design and specifications are subject to change without notice.
model number: TSX-304A MWPCT0008EAB 2023-11 CMSC/D/Printed in Japan

Canon Medical Systems Corporation meets internationally recognized
standards for Quality Management System ISO 9001, ISO 13485.
Canon Medical Systems Corporation meets the
Environmental Management System standard ISO 14001.

Catphan is a registered trademark of The Phantom Laboratory.

Made For life